



seca mBCA 525



Medical Body Composition Analyzer for determining body composition in a lying position

- + Specially designed and validated for mobile medical use
- + Patented measuring mat: BIA in supine position is finally comfortable, precise and easy
- + Lightweight and robust touchscreen monitor allows for direct evaluation
- + Unique accuracy through validation from a whole body MRI (skeletal muscle mass) and a 4C model (fat mass)
- + Validation studies are bundled and published in the European Journal of Clinical Nutrition 2013/2017, among others

mbca.seca.com

| | |
|-----------------------|--|
| Measurement method | 8-point Bioelectrical Impedance Analysis |
| Measurement current | 100 µA |
| Frequencies | 9 |
| Measured data storage | Up to 100,000 measurements |
| Measurement duration | 30 seconds |
| Software | seca analytics 115 |
| Accessories | + Rolling stand seca 475 + Basket seca 490 0019 for stand + Carrying case seca 432 |
| Consumables | + Single-use electrodes 490 0020 |

seca 959

EMR-validated chair scale with precise graduation

- + 50 g graduation for exact measurement while sitting
- + Fold-up arm and footrests for easy access
- + Sophisticated engineering and robust construction for a long lifespan
- + Suitable for patients affected by obesity with up to 300 kg weighing capacity



| | |
|----------------|--|
| Approval class | III |
| Capacity | 300 kg |
| Graduation | 50 g < 150 kg > 100 g < 250 kg > 200 g |
| Functions | TARE, Pre-TARE, HOLD, Auto-HOLD, BMI, Auto-CLEAR, auto weighing range switch-over, damping, automatic switch-off |
| Accessories | + Interface module seca 452 + Battery for longer mobile use seca 454 |

Medical Measuring Systems and Scales for

RENAL CARE



seca scales are indispensable tools for monitoring fluid balance, a critical aspect of treatment. By accurately measuring pre- and post-dialysis weight, seca scales enable healthcare providers to assess fluid removal efficacy, ensuring patients' safety and well-being throughout the process. seca offers a vast range of multifunctional scales to weigh immobile patients or patients with limited mobility as precise as possible.

seca 704

EMR-validated column scale with 300 kg capacity

- + Extremely robust, highly precise seca weighing technology with 300 kg capacity
- + Non-tilting base made of heavy duty cast iron
- + Impact- and scratch-resistant powder coating
- + Wide range of functions: BMI, pre-TARE, HOLD, Auto-HOLD, mother+child function, three-step damping
- + Expandable to measuring station with an optional height measuring rod

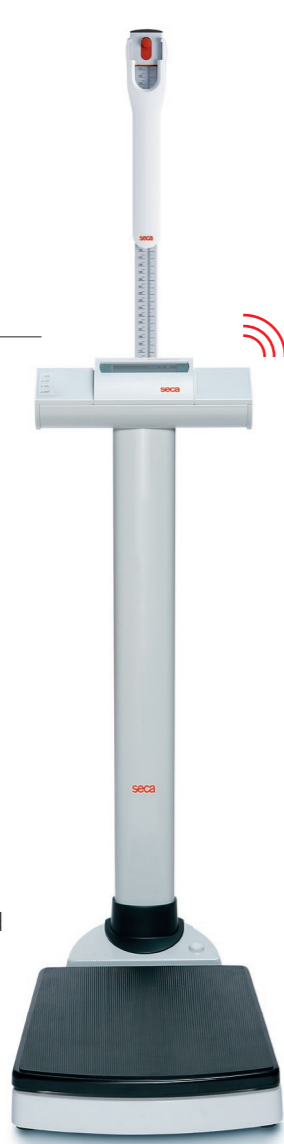


Image contains optional measuring rod seca 220.

Scale

| | |
|----------------|---|
| Approval class | III |
| Capacity | 300 kg |
| Graduation | 50 g < 150 kg > 100 g |
| Functions | TARE, Pre-TARE, mother+child function, HOLD, Auto-HOLD, BMI, Auto-CLEAR, auto weighing range switch-over, damping, automatic switch-off |
| Accessories | + Interface module seca 452 + Measuring rod seca 220 or panda seca 459 + Power adapter 447 + Switch-mode power adapter 400 |

Measuring rod seca 220 (optional)

| | |
|-----------------|-------------|
| Measuring range | 60 – 200 cm |
| Graduation | 1 mm |

seca 685

EMR-validated multifunctional scale

- + Safe weighing of standing, sitting or wheelchair patients
- + Large, easy to access platform
- + Wheel barrier provides extra safety
- + Integrated fold down seat
- + Easy to read display
- + Mobile due to integrated transport wheels



| | |
|----------------|--|
| Approval class | III |
| Capacity | 300 kg |
| Graduation | 100 g < 200 kg > 200 g |
| Functions | TARE, Pre-TARE, HOLD, Auto-HOLD, BMI, Auto-CLEAR, auto weighing range switch-over, damping |
| Accessories | + Interface module seca 452 |

seca 677

EMR-validated wheelchair scale with handrail

- + Large wheelchair accessible platform
- + Convenient and foldable to save space
- + Mobile due to transport wheels
- + Handrail as a standing aid
- + Non-slip rubber mat
- + Two ramps included



| | |
|----------------|--|
| Approval class | III |
| Capacity | 300 kg |
| Graduation | 100 g < 200 kg > 200 g |
| Functions | TARE, Pre-TARE, HOLD, Auto-HOLD, BMI, Auto-CLEAR, auto weighing range switch-over, damping |
| Accessories | + Interface module seca 452 |

seca 657

EMR-validated electronic platform scale with innovative memory function

- + Especially spacious platform for all common wheeled systems
- + Stable construction with very high capacity
- + Transportable due to wheels and integrated handles
- + Two ramps included
- + Innovative memory function for determining precise weight of patient



| | |
|----------------|--|
| Approval class | III |
| Capacity | 300 kg |
| Graduation | 100 g < 200 kg > 200 g |
| Functions | TARE, Pre-TARE, HOLD, Auto-HOLD, BMI, Auto-CLEAR, auto weighing range switch-over, damping, automatic switch-off |
| Accessories | + Interface module seca 452 + Stand for cable remote display seca 472 |

146.5 cm